

Sof'Airless



■ Airless piston technology



Nemera

Airless piston technology

SofAirless targets Rx and OTC dermal gels and creams and is already used for cold sore and acne treatments.

Can be used in any position thanks to the absence of air entry (360°)

For small volumes or reconstitution of unstable formulations



Upon pump actuation, the formulation is dispensed and the bottom piston rises. No air enters in contact with the formulation.

The **SofAirless** range is divided into two main families:

SofAirless XS



SofAirless XL



SofAirless XS

The solution for trial sizes or local treatments with very little skin surface to cover

Cap



SP24A pump



- Doses (μ l): 50, 80
- Fitment: Snap-on
- Actuator included




Piston & bottle



- Bottle volume: 2.05ml

SofAirless XS can be filled on a standard filling line. Our expert team is able to support you in implementing the automated process. We provide a manual snapping tool for initial trials.

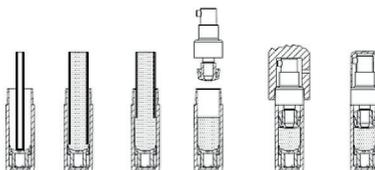

The pump design avoid air entrapment during snapping step.

SofAirless XL

The solution when drug reconstitution is needed

SofAirless XL* is delivered in 2 sub-assemblies

Dispensing head: SP943A+ & cap



- Dose: 1420 μ l
- Fitment: Screw
- Actuator included



Piston & bottle



- Bottle volume: 59ml

SofAirless XL can be filled on a standard filling line. The screwing fitment allows for formulation reconstitution at the pharmacy or at home by the patient.

* Not available for anti-acne applications

Dermal and transdermal delivery solutions for liquid and highly viscous formulations



Airless

Small Doses
(50µl - 200µl)

Sof'Airless XS

SP24A pump + piston bottle



Atmospheric

SP24, SP34

Pumps & actuators for viscous formulations:



SP27, SP270+, SP37, SP370+

Pumps & actuators for liquid formulations:

Sof'Bag®

SP943 pump + pouch-in-bottle sub-assembly



SP943

Pump & actuator for formulations with a wide range of viscosities:



Large Doses (672µl - 1550µl)

Sof'Airless XL

SP943A + pump + piston bottle



Delaminating systems

SP943A + delaminating bottles

Contact us for support on atmospheric bottle choice

Also available

Continuous valves for pressurized container CV20, 6668, 6685

Presented in this brochure

Contact : information@nemera.net
www.nemera.net

we put
patients
First